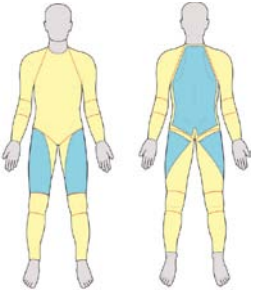
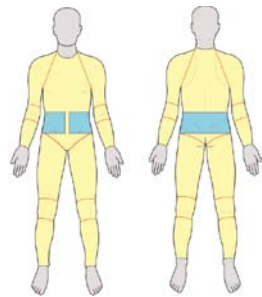
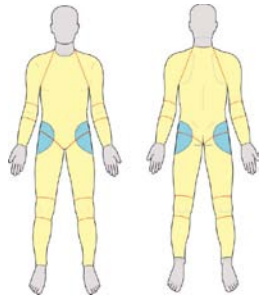
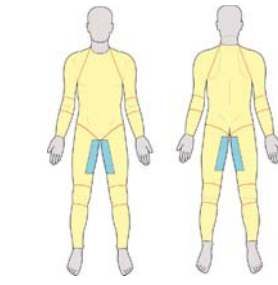
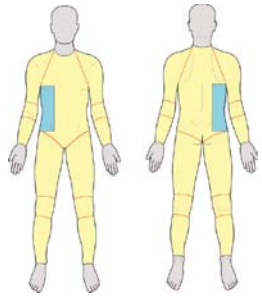
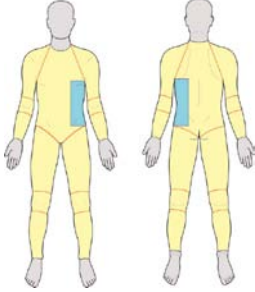
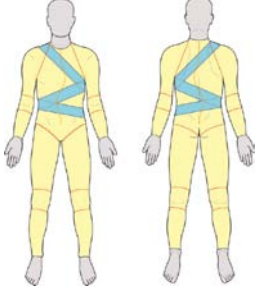
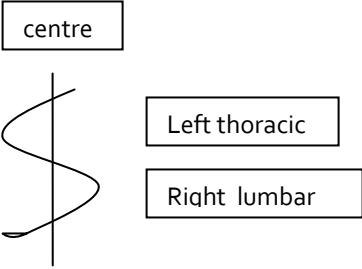
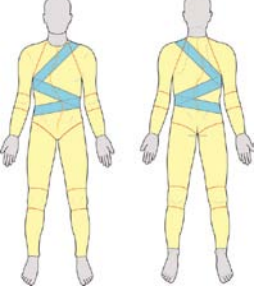
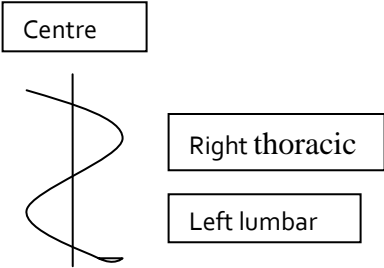










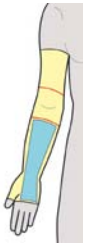

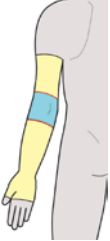






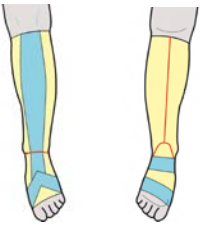
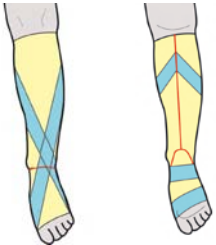
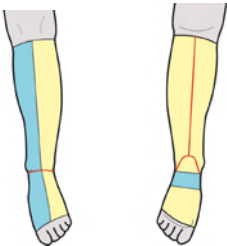
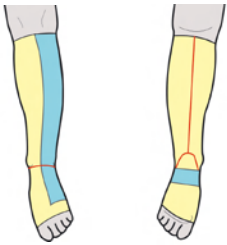
Code	Picture	Name and Description	Garment Type	Cost Category
SP		<p>Standard Panels</p> <p>To assist hip external rotation (turning the hips outwards) & back extension (straightening the back)</p>	<p>Suit Vest Leotard Pants Shorts</p>	Free of Charge
LP		<p>Lumbar Panel</p> <p>To assist postural control (ability to maintain upright posture)</p>	<p>Suit Vest Leotard</p>	Medium
HP		<p>Hip Panels</p> <p>To assist hip control (ability to tilt hips backwards and forwards) and external rotation (assist turning hips outwards)</p>	<p>Suit Pants Shorts</p>	Medium
ULM		<p>Medial Upper Leg Panels</p> <p>To assist adduction (bringing legs together) and internal rotation (assist turning hips inwards)</p>	<p>Suit Pants Shorts</p>	Medium
TPL		<p>Lateral Trunk Panels</p> <p>To assist correction of left lateral flexion (panels positioned on right to pull the trunk from the left)</p>	<p>Suit Leotard</p>	Medium

<p>TPR</p>		<p>Lateral Trunk Panels</p> <p>To assist correction of right lateral flexion (panels positioned on left to pull the trunk from the right)</p>	<p>Suit Leotard</p>	<p>Medium</p>
<p>SPR</p>		<p>Scoliosis Panels</p> <p>To assist correction of left thoracic (area of upper back) and right lumbar (area of lower back) curves</p>  <p>centre</p> <p>Left thoracic</p> <p>Right lumbar</p>	<p>Suit</p>	<p>Complex</p>
<p>SPL</p>		<p>Scoliosis Panels</p> <p>To assist correction of right thoracic (area of upper back) and left lumbar (area of lower back) curves</p>  <p>Centre</p> <p>Right thoracic</p> <p>Left lumbar</p>	<p>Suit</p>	<p>Complex</p>
<p>AP</p>		<p>Anterior Abdominal Panels</p> <p>To assist trunk flexion (assist the trunk to bend forwards)</p>	<p>Suit Leotard</p>	<p>Medium</p>

BP		<p>Posterior Back Panels</p> <p>To assist trunk extension (assist the trunk to straighten)</p>	<p>Suit Leotard</p>	<p>Medium</p>
KE		<p>Anterior Leg Panels</p> <p>To assist knee extension (to assist straightening of the knee)</p>	<p>Suit</p>	<p>Complex</p>
KF		<p>Posterior Leg Panels</p> <p>To assist knee flexion (to assist bending of the knee)</p>	<p>Suit</p>	<p>Complex</p>
APP		<p>Anterior Pelvic Panels</p> <p>To assist posterior tilt (to assist the pelvis to rock backwards)</p>	<p>Suit Pants Shorts</p>	<p>Complex</p>
PPP		<p>Posterior Pelvic Panels</p> <p>To assist anterior tilt (to assist the pelvis to rock forwards)</p>	<p>Suit Pants Shorts</p>	<p>Complex</p>

<p>CP1</p>		<p>Anterior Chest Panels</p> <p>To assist protraction of the scapulae (to assist rounding forwards of the shoulders)</p>	<p>Suit Vest Leotard</p>	<p>Medium</p>
<p>CP2</p>		<p>Posterior Cross Panels</p> <p>To assist retraction of the scapulae and back extension (to assist taking the shoulders back and straightening of the back)</p>	<p>Suit Leotard Vest</p>	<p>Medium</p>
<p>HE</p>		<p>Posterior Leg Panels</p> <p>To assist hip extension (to assist straightening of the hips)</p>	<p>Suit</p>	<p>Complex</p>
<p>UL1</p>		<p>Dorsal Wrist Panel</p> <p>To resist (to assist in preventing the wrist being bent forwards) wrist flexion (to assist wrist extension)</p>	<p>Sleeve Gauntlet</p>	<p>Free of Charge</p>
<p>UL2</p>		<p>Volar Wrist Panel</p> <p>To resist (to assist in preventing the wrist being bent backwards) wrist extension (to assist wrist flexion)</p>	<p>Sleeve Gauntlet</p>	<p>Minor</p>

UL3		<p>Anterior Elbow Panel</p> <p>To resist elbow hyperextension (elbow being bent backwards beyond straight)</p>	Sleeve	Minor
UL4		<p>Posterior Elbow Panel</p> <p>To resist elbow flexion (elbow being bent forwards)</p>	Sleeve	Minor
UL5		<p>Posterior Spiral Forearm Panel BACK</p> <p>To resist pronation (to assist the forearm to be turned palm upwards)</p>	Sleeve Gauntlet	Complex
UL6		<p>Anterior Spiral Panel FRONT</p> <p>To resist supination (to assist the forearm to be turned palm downwards)</p>	Sleeve Gauntlet	Complex
UL7		<p>Thumb Panel</p> <p>To assist extension & abduction (to assist thumb straightening and taking out of the palm)</p>	Sleeve Gauntlet	Minor
UL8		<p>Diagonal Wrist Panel</p> <p>To correct ulnar deviation (to stop the wrist from falling sideways away from the thumb)</p>	Sleeve Gauntlet	Minor

<p>LL1</p>		<p>Lower Leg Central Panel</p> <p>To assist dorsi-flexion (to assist lifting the foot upwards) and control alignment (keep straight)</p>	<p>Sock</p>	<p>Complex</p>
<p>LL2</p>		<p>Lower Leg Cross Panel</p> <p>To assist dorsi-flexion (to assist lifting the foot upwards)</p>	<p>Sock</p>	<p>Complex</p>
<p>LL3</p>		<p>Lower Leg Medial Panel</p> <p>To assist dorsi-flexion (to assist lifting the foot upwards) and resist eversion(pronation) (to assist in stopping the foot turning down and outwards)</p>	<p>Sock</p>	<p>Medium</p>
<p>LL4</p>		<p>Lower Leg Lateral Panel</p> <p>To assist dorsi-flexion (to assist lifting the foot upwards) and resist inversion(supination) (to assist in stopping the foot turning down and inwards)</p>	<p>Sock</p>	<p>Medium</p>